



# The NEW Micro+™

Helping people to stop smoking, one breath at a time.



CE  
0086




One CO monitor for all your testing needs

**FREE**  
COdata+™  
software

Ideal for in-depth stop smoking advice delivered by:

- Stop smoking services
- General practice
- Schools
- Pharmacies
- Health visitors
- Midwives
- Researchers

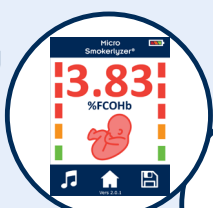

Automatic %FCOHb conversion



For ambient levels of CO



Create and save patient details

**Smokerlyzer®**  
[www.bedfont.com](http://www.bedfont.com)

# Features and Benefits

- Brand new easy to use interface
- Create and store patient details
- SteriTouch® technology to kill off harmful bacteria
- Large touch screen
- Automatic conversion to %FCOHb
- Familiar green, amber and red traffic light system
- Ambient levels of CO

## Technical Specification

<b>Concentration Range (CO)</b>	0-500ppm
<b>Display</b>	Full colour touchscreen
<b>Detection Principle</b>	Electrochemical sensor
<b>Repeatability</b>	<3% difference on consecutive readings
<b>H<sub>2</sub> Cross Interference</b>	<6%
<b>Power</b>	3 x AA (LR6 or Equivalent) 1 x CR2032 Lithium coin cell
<b>T90 Response Time</b>	≤17 seconds
<b>Operating Temperature Range</b>	10-40°C
<b>Operating/Transport/Storage Pressure</b>	800-1200mbar
<b>Operating Humidity</b>	10-90% non-condensing
<b>Sensor Operating Life</b>	5 years – 5 year warranty
<b>Sensor Sensitivity</b>	1ppm
<b>Sensor Drift</b>	<10% per annum
<b>Dimensions</b>	Approx. 34 x 75 x 140 mm
<b>Weight</b>	Approx. 200g (including batteries)
<b>Monitor Construction</b>	Case: Polycarbonate/ABS blend SteriTouch anti-Microbial additive D-Piece: Polypropylene



## British Technology

### Bedfont Scientific Ltd

Station Road, Harrietsham, Maidstone  
Kent, ME17 1JA, England  
Tel: +44 (0)1622 851122, Fax: +44 (0)1622 854860  
Email: [ask@bedfont.com](mailto:ask@bedfont.com)  
[www.bedfont.com](http://www.bedfont.com)

© Bedfont Scientific Limited 2015

Issue 5- March 2016- Part No: MKT375  
Bedfont Scientific Limited reserve the right to change or update this literature without prior notice.  
Registered in: England and Wales. Registered No: 1289798

